EN

MASSFORCE

MASSEYES

model 10S-C

INSTRUCTION MANUAL

Thanks for your purchase.

For any inquiries, contact us at hello@massforce.co or on our website www.massforce.co

Read Instructions
Before Operating

Imported by OUI SMART
15 RUE DES HALLES
75001 PARIS, FRANCE
www.massforce.co

FC

CE

 **RoHS**

WARNINGS

- This product has a heating function. For users who are prone to heat allergies or insensitive to heat sources, please use the heating function with caution.
- If you feel unwell during or after use, please stop using this product.
- This product is a personal care product and cannot replace professional treatment or health care. Please do not use it for diagnosis or other commercial purposes.
- Do not use it with other medical equipment or more than 2 sets of this product at the same time.
- Do not use the same site for more than 30 minutes, otherwise, it may cause discomfort.
- Do not use the heating function for more than 30 minutes, otherwise, it may cause low-temperature skin burns.
- Do not use the product after accidentally falling into water or liquid.

PRECAUTIONS

- Using time of 15 minutes is recommended. Overtime use may cause discomfort.
- Please do not use the product immediately after use eye drops.
- Please do not use with facial mask together.
- If you feel over hot, please turn OFF the product.
- This product contains lithium batteries, please do not disassemble, modify, immerse the product in water or throw fire into it by yourself.
- In the event of an abnormal situation that cannot work, please cut off the power immediately and remove it, check the basic troubleshooting or consult customer service.
- When the product is working, do not touch the electrode pads with metal accessories.
- Do not charge, use, and place in high-temperature places such as near the stove or under the scorching sun.
- This product is not waterproof, please do not wash and pay attention to waterproof.
- This product is can only be decomposed when the product is discarded. If you need to repair it, please consult customer service in time.
- Do not leave the product in front of unattended children.
- Do not allow children or persons with limited ability to use this product is not suitable for any sick or incapacitated person.
- If you have one of the following conditions, please consult a doctor before use:
 - A. Those are pregnant, sick or in poor health,
 - B. Patients who are undergoing treatment and feel abnormalities,
 - C. Patients with high fever and infectious diseases,
 - D. Patients with hypertension, heart disease, cerebral thrombosis, arteriosclerosis or malignant tumors,
 - E. People with skin sensation disorder or abnormal skin sensitivity,
 - F. Those in the postoperative recovery stage,
 - G. Those who are uncertain about their health.

ADVERSE EFFECTS

- On very rare occasions, first-time users of EMS have reported feeling light-headed or faint. We recommend that you use the product while seated until you become accustomed to the sensation.

You should stop using the device and consult with your physician if you experience adverse reactions from the device.

SPECIFICATION

Product model: 10S-C
Power input: DC 5V (USB charging port)
Battery: 850 mAh
Charging time: 2 hours
Working time: 60-90 minutes

Package dimensions: 19 x 14 x 8 cm;
Weight 0.5 kg
Color: white
Heating function: 38°C-42°C

What's in the box:
1 unit Masseyes
1 unit USB cable
1 pcs user manual

Easy button control
User friendly

Ergonomic curving
It doesn't press on the nose

Adjustable tightness
You can adjust the comfortable size

Soft inside
Comfortable and skin-friendly

***Note:**

This product is not suitable for medical equipment.

Does not have a therapeutic effect, children can not touch.

INTRODUCTION

Masseyes is utilizing advanced micro-controlled steaming technology for eye reflexology, relief of dry eyes, tired eyes, and fatigue. This adopts its thermostat-regulated soothing heat wave with its warm steam to penetrate the skin tissues around the eyes, and thus the thermal reaction between the molecules accelerates blood circulation and helps relieve eye fatigue. Its built-in rechargeable battery and touch-panel control let you sit back with peace of mind and enjoy its luxury relief of dry eyes, tired eyes, and fatigue.

Well-designed for myopia, long-sightedness, astigmatism, eye fatigue. This is good for people who have fine lines around the eyes, and dark circles. People who have eye ache, headache, psychasthenia, insomnia.

180° foldable design allows you to pack it into your bag and make a quick massage possible anywhere with adjustable elastic straps bring you proper size. Masseyes is built-in heating pad that offering a moderate temperature between 38°C-42°C, and heating up within 1 or 2 seconds to stimulate your blood circulation around eyes.


USING MASSEYES™

Charging mode:

- Fold the product first, and use the USB cable to charge
- Please use 5V/1A-2A regular adapter to charge the product. DO NOT USE ANY ADAPTER WITH OVER 5V
- Please check the indicator light. Please do not use the product while it is being charged. When the product makes a "Du Du" sound, please charge the machine. Long-time use of the product with a low battery may cause a charging problem.

Host

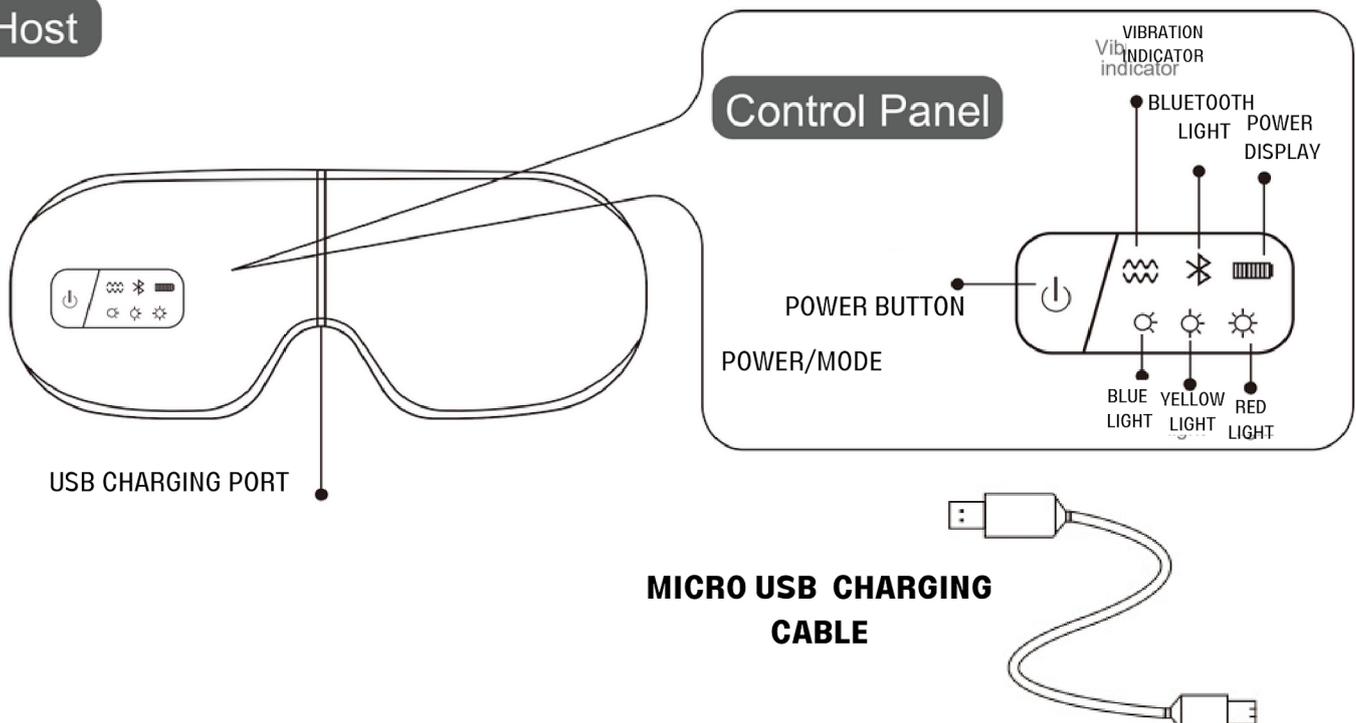**ON/OFF:**

Press the power button for 3 seconds to turn ON the Masseyes
To turn it OFF, long press for 3 seconds.

Masseyes heats up in seconds, providing immediate warmth and soothing relief. It's perfect for relieving eye strain, headaches, and migraines, sinus pressure, post-nasal drip, dryness and more!

**Thermotherapy function:**

Low-temperature warm compress, Blue light display

- Press the power button one more time, it will turn to Medium -temperature warm compress, Yellow light display
- Press the power button one more time, it will turn to High-temperature warm compress, Red light display
- Press the power button one more time, it will turn to the Highest temperature warm compress

his eye massager comes with vibration mode to help relieve stress and tension around the eyes. The gentle massage will help increase blood circulation and reduce puffiness and dark circles.

Vibration function:

1. Double-click to change to vibration 1 for acupressure mode
2. Click again to convert to vibration 2 for eye relaxation mode
3. Double-click again to return to the original mode

MAINTENANCE

Please do not use it when charging, start Masseyes after you wear it perfectly to prevent excessive pressure on the airbag, the damage it.

Do not soap or rinse in water. Water or other liquids should be prevented from entering the body during the cleaning process to avoid short circuit and electric shock hazards during use.

Use a dry and soft cloth or paper to clean the product.

Do not clean it while it is charging.

TROUBLESHOOTING

If you encounter problems when using the product, please refer to the following tips and solutions. If the product is still in a faulty state or cannot be used, please turn off the power.

Users are not allowed to disassemble or repair the product by themselves. The product damage caused by disassembly or repair by themselves will invalidate the product warranty, and the company will not bear any consequences caused by this situation.

Unable to start - The internal charging battery is low. Please confirm that the key of the host product is pressed.

No charging or all the indicator light, not shiny when charging - the product is in the low battery due to did not use for a long time, full charge the product.

All the indicator lights shiny or There is "Di Di" sound appears - the product is in the low battery due to did not use for a long time, full charge the product.

There is a Siren when charging - the product is in the low battery due to did not use for a long time, full charge the product.

When you feel insufficient temperature when under warm mode - low battery, full charge the product.

When air pressure mode and warm mode unable to start - Caused by the status of power supply, full charge the product.

Barometric intensity too strong or weak - Adjust the headband to the length you need.

Murmur sound appears when under air pressure mode - When the airbag start, there would be murmur sound appear, no need to take any action.

How to judge charging situation - The indicator light keeps shiny when charging. When the product full charged, the indicator turns off.

Massforce branded items purchased from massforce.co include a 2-year limited replacement warranty. For Customer Service, please contact us via email hello@massforce.co.

DO NOT DISPOSE OF BATTERIES IN FIRE. ALWAYS DISPOSE OF BATTERIES PER LOCAL & FEDERAL GUIDELINES. The Massforce® name and logo are registered trademarks.

©Massforce. All Rights Reserved.

Read Instructions
Before Operating

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