
EN

MASSFORCE

MASSFIT™ Abdominal Belt

model MH1082

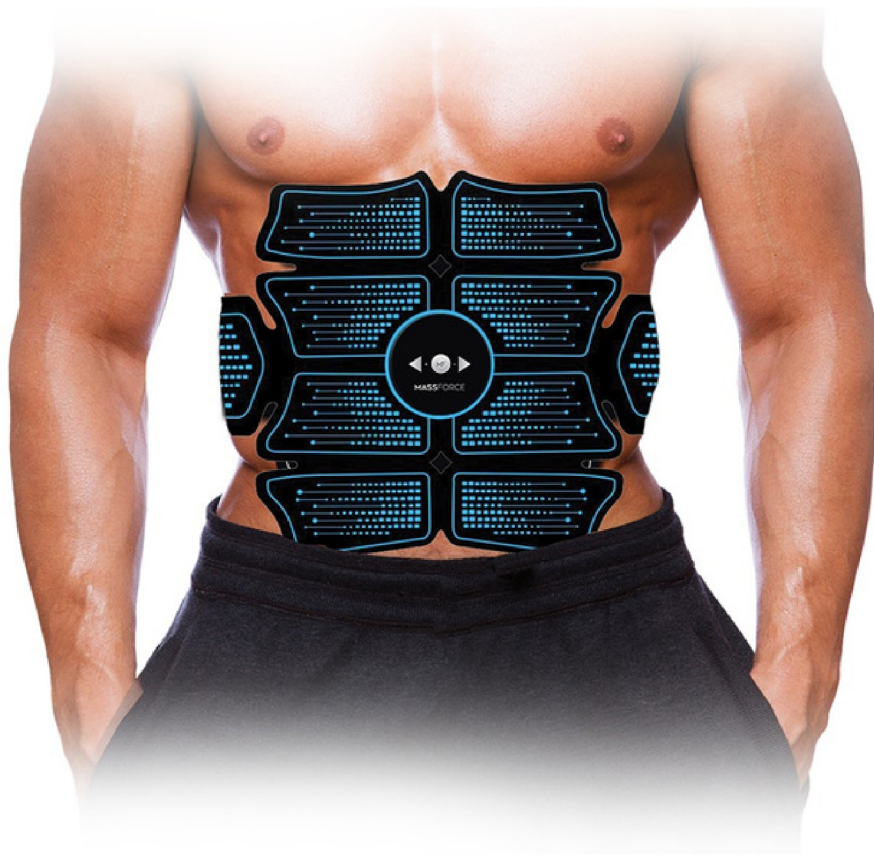INSTRUCTION MANUAL

Read the manual before operating this product.

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www.massforce.co

WARNINGS

- If you are in the care of a physician, consult your physician before using this device.
- DO NOT use this device if you have a cardiac pacemaker, implanted defibrillator, or other implanted metallic or electronic devices.
- DO NOT place this device over your neck.
- DO NOT place this device across your chest.
- DO NOT use this device during pregnancy.
- DO NOT use this device on children.
- DO NOT place this device over the carotid sinus nerves, the front of the neck, or around the mouth.
- DO NOT open unit. The opening will void the warranty.
- DO NOT use this device over, or in proximity to, cancerous lesions.
- DO NOT place the device over open wounds, sores, or rashes, or over swollen, red, infected, or inflamed areas or skin eruptions (e.g., phlebitis, thrombosis, varicose veins).
- DO NOT use this device when in bath or shower.
- DO NOT use this device while sleeping.
- DO NOT use this device while driving or operating machinery.
- DO NOT use this device in the presence of electromagnetic fields.
- DO NOT use this device if unit has been damaged.
- DO NOT use on children or pets.
- DO NOT use in an area where water or moisture can get on the unit, such as a bathroom or other places with high humidity.

PRECAUTIONS

- The long-term effects of electrical stimulation are unknown.
- Since the effects of stimulation of the brain are unknown, stimulation should not be applied across your head.
- The safety of electrical stimulation during pregnancy has not been established.
- Use caution if you have a tendency to bleed internally, such as following an injury or fracture.
- You may experience skin irritation or hypersensitivity.
- If you have suspected or diagnosed heart disease, you should follow the precautions recommended by your physician.
- Keep this device out of the reach of children.
- If you have suspected or diagnosed epilepsy, you should follow the precautions recommended by your physician.
- Consult with your physician prior to using the device after a recent surgical procedure, because the stimulation may disrupt the healing process.
- Use caution if stimulation is applied over the menstruating or pregnant uterus.
- Use caution if stimulation is applied over areas of skin that lack normal sensation.
- Use this device only with the electrodes and accessories recommended by the manufacturer.
- Avoid placing the pads over metal implants.

ADVERSE EFFECTS

- Skin irritation and burns beneath the stimulation electrodes applied to your skin can occur.
- On very rare occasions, first-time users of EMS have reported feeling light-headed or faint. We recommend that you use the product while seated until you become accustomed to the sensation.

You should stop using the device and consult with your physician if you experience adverse reactions from the device.

SPECIFICATION

Power Supply: DC 3.7v
Output Voltages: 55 V at 500Ω
Pulse Width: 100 μS
Frequency: 0-200 Hz
Timer: 20 min
Intensity: 10 levels
Modes: 6 modes

Operation Condition: -10~40 °C, 30%~85% Humidity
Storage Condition: -10~50 °C, 10%~95% Humidity
Transportation Condition: -10~50 °C, 35%~85% Humidity


INTRODUCTION

MASSFIT™ Abdominal Belt delivers electrical pulses through the included electrodes. Each features multiple modes of different pulse frequencies, providing Electrical Muscle Stimulation (EMS), as well as Transcutaneous Electrical Nerve Stimulation (TENS). There are 6 pulse modes, each with 10 intensity levels. Each unit includes conductive gel pads and one USB charging cable.

Indications for Use

Use Abdominal Belt to stimulate healthy abdominal muscles in order to improve and facilitate muscle performance. It is also intended for temporary relief of pain associated with sore and aching abdominal muscles from exercise or normal work activities.

UNIT DIAGRAM

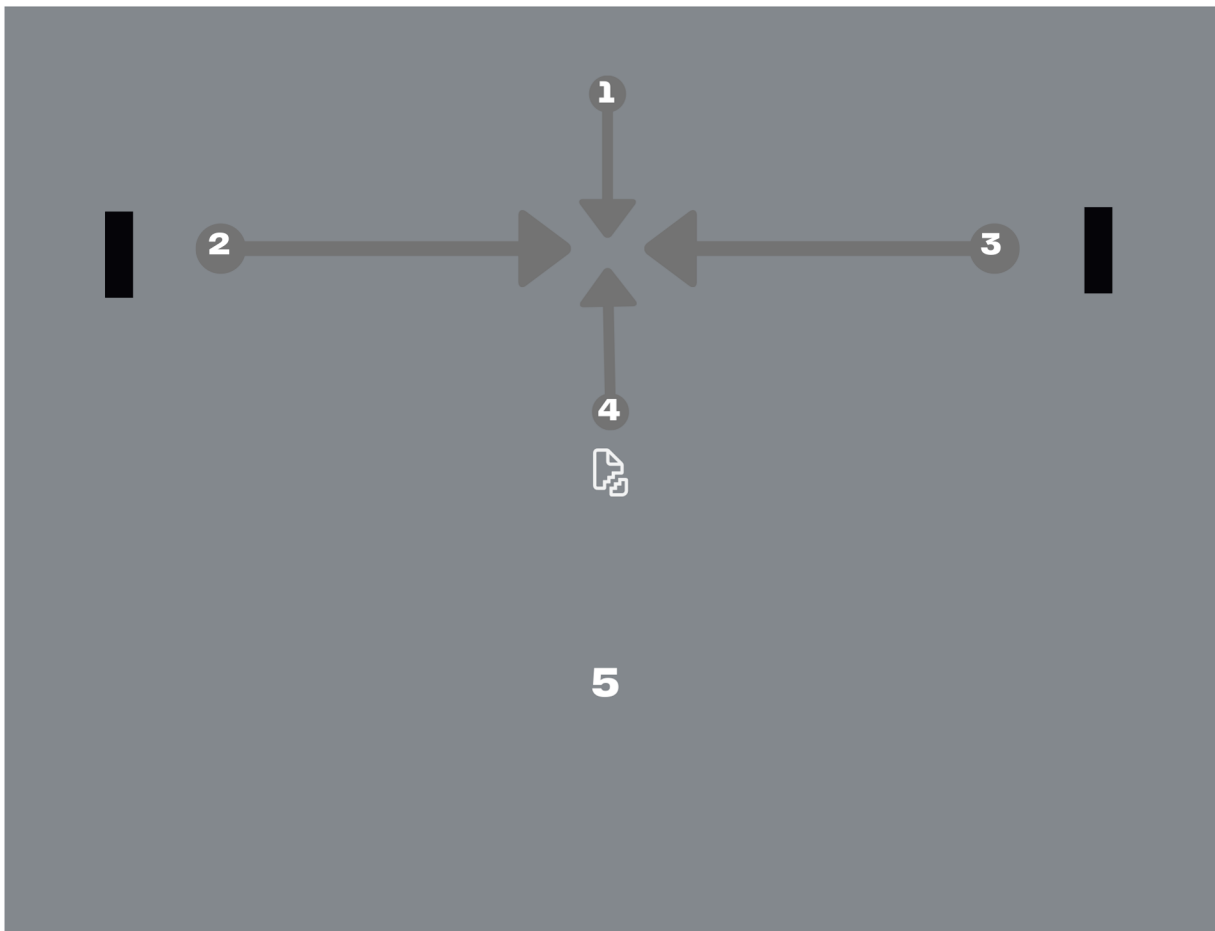1. Adjustment button "M". You can choose a mode and a different intensity to meet your needs.
2. Plus the "+" button to increase the intensity.
3. Minus the "-" button to increase the intensity. 2x quick press for turning OFF.
4. Light indicator to show you which mode is ON.
5. Electrodes on gel pads (on the back of the unit).

USING MASSFIT™

Abdominal Belt

1. Place the device to MASSFIT™ Abdominal Belt by clicking two pins.
2. Making sure your skin is clean and dry, place the MASSFIT™ Abdominal Belt onto the muscle group you wish to target with the MASSFIT™ Abdominal Belt turned off.
3. If using the Massforce Abdominal Belt, the device should be located over the belly button.
4. Turn the electrodes on by pressing the “M” power button. The intensity will be in the lowest setting. There are 6 modes you can adjust and easily be seen by the light indicator.
5. Press the + or - buttons to adjust intensity. There are 10 intensity levels.
6. When the session is complete, power off the unit by pressing the Minus “-“ button with a 2x quick press.
7. To charge the unit, connect the USB cable to the device and plug it into a USB wall adapter for 20 minutes.
8. Best results are achieved with a combination of diet and exercise.

Use for 20 minutes at a time, 1-2 times per day per area.

PULSE MODES

The device has multiple modes, covering EMS, as well as TENS. If you are using the device for the first time, it is recommended that you start with the default Mode 1 for improving muscle performance. After the muscle training, you may feel sore and aching muscles due to strain from exercise. We recommend you follow up with one of the TENS modes, as they are good for the relief of such pain. Please see the EMS and TENS modes in the following table.

10 Intensities						
	LV 1		LV 5		LV 10	
6 Mode						
	A	B	C	D	E	F

Easily click the “M” button to switch the mode and “+” to increase the intensity level. The mode will change be indicated by the LED light flashing

1. Mode 1 (EMS): Stimulation of healthy muscles in order to improve and facilitate muscle performance
2. Mode 2 (TENS): Relief of pain associated with sore and aching muscles due to strain from exercise or work activities
3. Mode 3 (EMS and TENS): Combination of the benefit from both EMS and TENS of Mode 1 & 2
4. Mode 4 (EMS): Visible muscle contraction saw for improving muscle performance
5. Mode 5 (TENS): Deeper relief of pain associated with sore and aching muscles due to strain from exercise or work activities
6. Mode 6 (EMS and TENS): Combination of the benefit from both EMS and TENS of Mode 4 & 5

Modes and functions:

- Mode A & B: Fat burning.

In these two modes, the muscles move intensely and burn fat quickly.

- Mode C & D: muscle growth.

Strong elastic movement, deep stimulation to support muscle growth.

- Mode E & F: body sculpting.

These two modes are composite modes that help build a perfect body shape.

GEL PAD MAINTENANCE

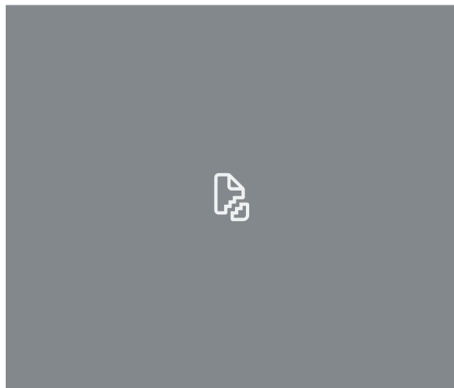
Figure 1

- Before applying the device to the skin, make sure to clean the desired area beforehand. To cleanse the skin, simply wet a towel with mild soap and water and gently wipe the area. This helps to remove any natural oils that can remove the adhesion of the electrodes.

- It is recommended to remove excess hair from the desired training area in order to optimize the life of the product.

- Do not apply lotions or oils before using the device, as this will prevent the gel pad from adhering to the skin.

- When you are finished with your workout, a cleanse is in order. To do this, use a damp cloth or towel and clean them gently. Then be sure to place the device on the protective plastic sheet.

- Use only the plastic sheet provided to store the product and gel pads. This sheet is designed to maximize the life of gel pads.

- The lifespan of gel cushions varies depending on the frequency of use and the care given. To extend their life, keep the gel pads clean and store them properly on the plastic sheet provided.


- After several uses, the contact gel on your electrodes may dry out. Sometimes it is enough to moisten the electrodes to reactivate this gel by adding a few drops of water or by wiping lightly with a damp towel. Be careful, this gesture should not be repeated too frequently, otherwise it will have the opposite effect and the electrodes will lose their adhesion.

- Replace the gel pads when they lose their stickiness.

To remove the electrodes from the gel pads, gently roll the edge of the gel pad with your thumb. If you roll it properly, the gel pad will come off in one piece, as shown in figure 1.

TROUBLESHOOTING

- Stimulation is weak or non-existent: Be sure skin is clean and pads are firmly attached.
- Device does not turn on: Recharge for 20 minutes.
- Gel pads are not sticky: The gel pads will lose adhesiveness over the course of use. Please visit [Massforce.co](https://massforce.co) for replacements.

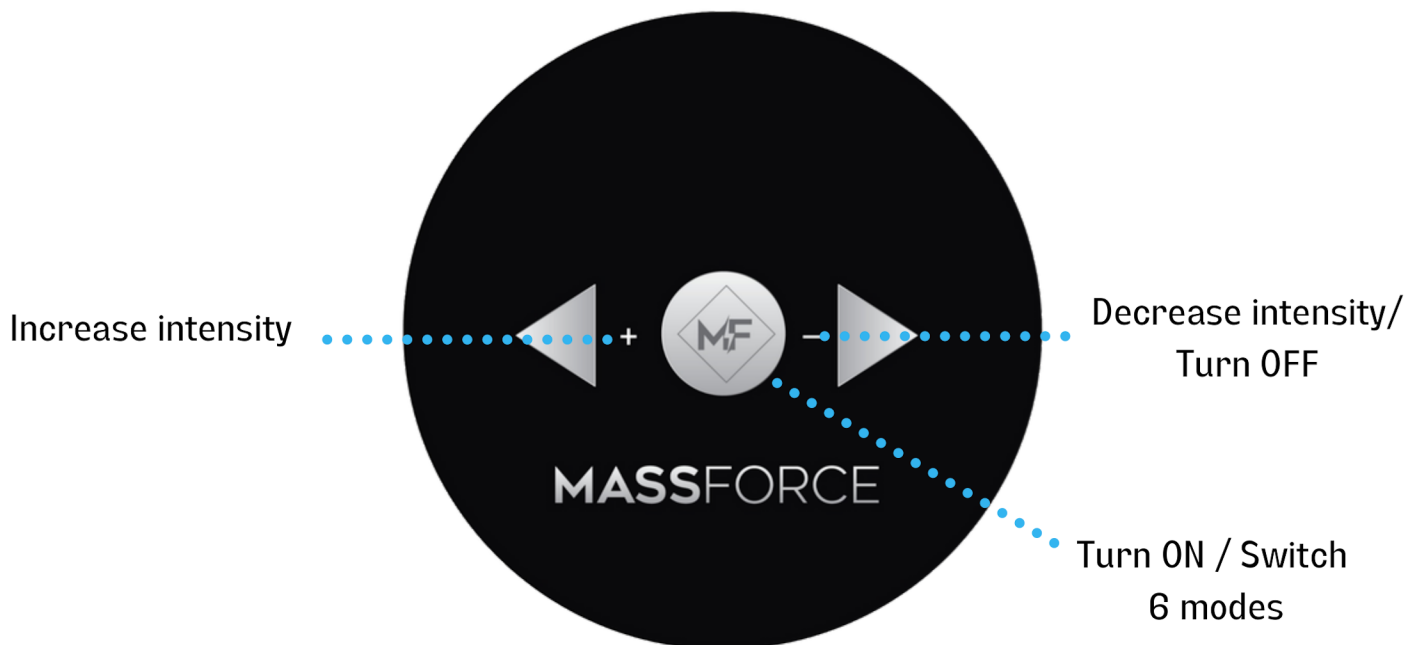Massforce branded items purchased from massforce.co include a 1-year limited replacement warranty. For Customer Service, please send us email: Hello@massforce.co.

DO NOT DISPOSE OF BATTERIES IN FIRE. ALWAYS DISPOSE OF BATTERIES PER LOCAL & FEDERAL GUIDELINES. The SHARPER IMAGE® name and logo are registered trademarks.


©Massforce. All Rights Reserved.

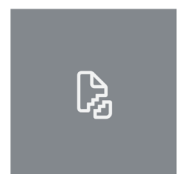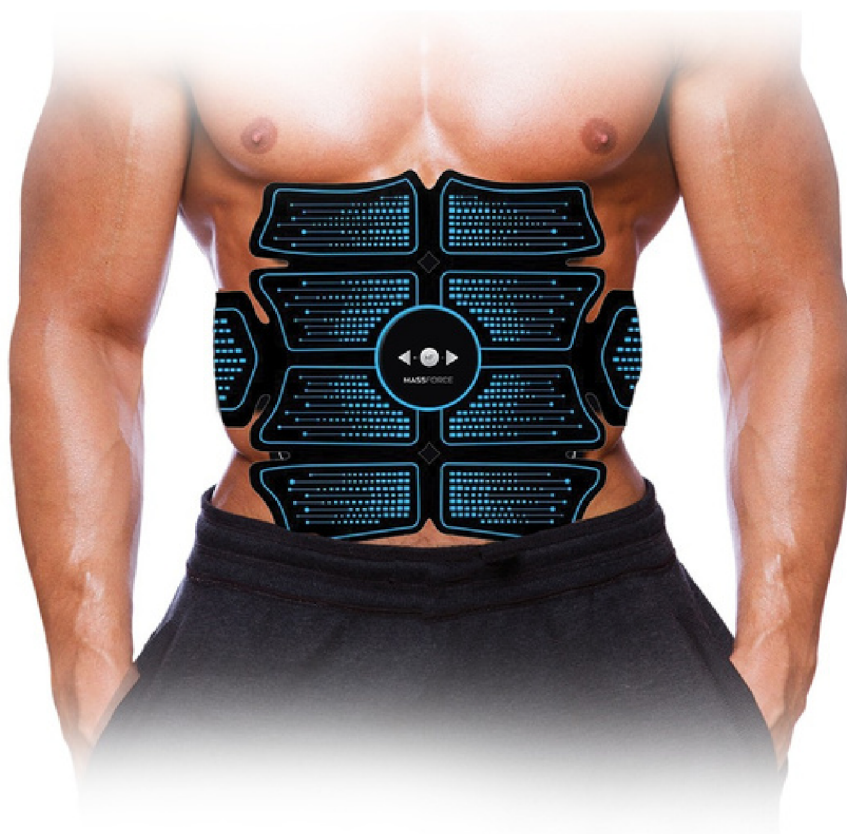
Read Instructions
Before Operating

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