FR

MASSFORCE

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Model MH-2017

Thanks for your purchase.

For any inquiries, contact us at hello@massforce.co or on our website www.massforce.co

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15 RUE DES HALLES
75001 PARIS, FRANCE

www.massforce.co

Read Instructions
Before Operating

CE

Made in China

WARNINGS

- If you are in the care of a physician, consult your physician before using this device.
- DO NOT use this device if you have a cardiac pacemaker, implanted defibrillator, or other implanted metallic or electronic devices.
- DO NOT place this device over your neck.
- DO NOT place this device across your chest.
- DO NOT use this device during pregnancy.
- DO NOT use this device on children.
- DO NOT place this device over the carotid sinus nerves, the front of the neck, or around the mouth.
- DO NOT open unit. The opening will void the warranty.
- DO NOT use this device over, or in proximity to, cancerous lesions.
- DO NOT place the device over open wounds, sores, or rashes, or over swollen, red, infected, or inflamed areas or skin eruptions (e.g., phlebitis, thrombosis, varicose veins).
- DO NOT use this device when in bath or shower.
- DO NOT use this device while sleeping.
- DO NOT use this device while driving or operating machinery.
- DO NOT use this device in the presence of electromagnetic fields.
- DO NOT use this device if unit has been damaged.
- DO NOT use on children or pets.
- DO NOT use in an area where water or moisture can get on the unit, such as a bathroom or other places with high humidity.

PRECAUTIONS

- The long-term effects of electrical stimulation are unknown.
- Since the effects of stimulation of the brain are unknown, stimulation should not be applied across your head.
- The safety of electrical stimulation during pregnancy has not been established.
- Use caution if you have a tendency to bleed internally, such as following an injury or fracture.
- You may experience skin irritation or hypersensitivity.
- If you have suspected or diagnosed heart disease, you should follow the precautions recommended by your physician.
- Keep this device out of the reach of children.
- If you have suspected or diagnosed epilepsy, you should follow the precautions recommended by your physician.
- Consult with your physician prior to using the device after a recent surgical procedure, because the stimulation may disrupt the healing process.
- Use caution if stimulation is applied over the menstruating or pregnant uterus.
- Use caution if stimulation is applied over areas of skin that lack normal sensation.
- Use this device only with the electrodes and accessories recommended by the manufacturer.
- Avoid placing the mat over metal implants.


ADVERSE EFFECTS

- Skin irritation and burns beneath the stimulation electrodes applied to your skin can occur.
- On very rare occasions, first-time users of EMS have reported feeling light-headed or faint. We recommend that you use the product while seated until you become accustomed to the sensation.

You should stop using the device and consult with your physician if you experience adverse reactions from the device.

SPECIFICATION

Power Supply:	DC 3.7v	
Output Voltages:	55 V at 500Ω	
Pulse Width:	100 μS	Operation Condition: -10~40 °C, 30%~85% Humidity
Frequency:	0-200 Hz	Storage Condition: -10~50 °C, 10%~95% Humidity
Timer:	20 min	Transportation Condition: -10~50 °C, 35%~85% Humidity
Intensity:	10 levels	
Modes:	6 modes	

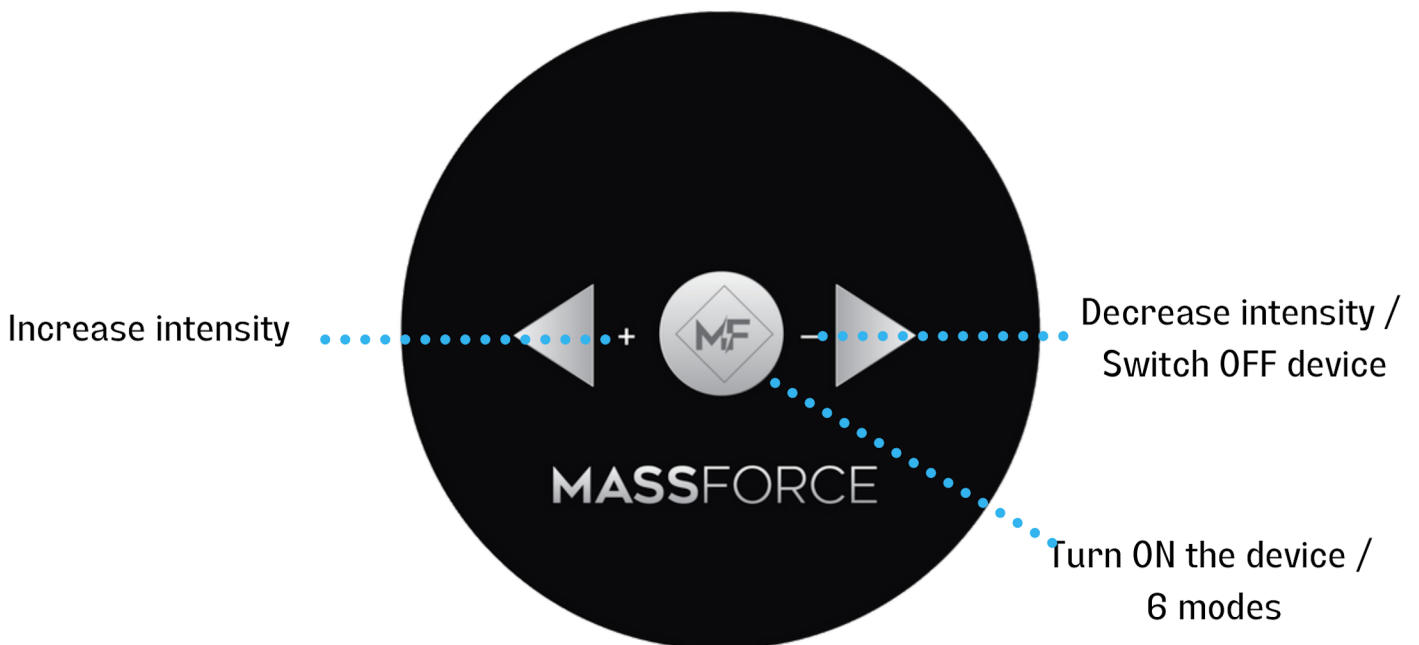INTRODUCTION

MASSFOOT™ EMS Tens Foot Massager delivers electrical pulses through the included electrodes. Each features multiple modes of different pulse frequencies, providing Electrical Muscle Stimulation (EMS), as well as Transcutaneous Electrical Nerve Stimulation (TENS). There are 6 pulse modes, each with 10 intensity levels. Each unit includes conductive mat and one USB charging cable.

Indications for Use

Use MASSFOOT™ EMS Tens Foot Massager to stimulate healthy leg muscles in order to improve and facilitate muscle performance. It is also intended for temporary relief of pain associated with sore and aching leg muscles from exercise or normal work activities.

UNIT DIAGRAM

1. Adjustment button "M". You can choose a mode and a different intensity to meet your needs.
2. Plus the "+" button to increase the intensity.
3. Minus the "-" button to increase the intensity. 2x quick press for turning OFF.
4. Light indicator to show you which mode is ON.

USING MASSFOOT™

EMS Tens Foot Massager

1. Place the device to MASSFOOT™ EMS Tens Foot Massager by clicking two pins.
2. Making sure your sole of the feet are clean and dry, place the MASSFOOT™ EMS Tens Foot Massager on the floor turned off.
3. If using the MASSFOOT™ EMS Tens Foot Massager, the device should be located on a flat surface.
4. Turn the electrodes on by pressing the “M” power button. The intensity will be in the lowest setting. There are 6 modes you can adjust and easily be seen by the light indicator.
5. Press the + or - buttons to adjust intensity. There are 10 intensity levels.
6. When the session is complete, power off the unit by pressing the Minus “-” button with 2x quick press.
7. To charge the unit, connect the USB cable to the device and plug it into a USB wall adapter for 20 minutes.
8. Best results are achieved with a combination of diet and exercise.

Use for 20 minutes at a time, 1-2 times per day per area.

PULSE MODES

The device has multiple modes, covering EMS, as well as TENS. If you are using the device for the first time, it is recommended that you start with the default Mode 1 for improving muscle performance. After the muscle training, you may feel sore and aching muscles due to strain from exercise. We recommend you follow up with one of the TENS modes, as they are good for the relief of such pain. Please see the EMS and TENS modes in the following table.

10 Intensities										
	LV 1		LV 5						LV 10	
6 Mode										
	A	B	C	D	E	F				

Easily click the “M” button to switch the mode and “+” to increase the intensity level. The mode will change be indicated by the LED light flashing

1. Mode 1 (EMS): Stimulation of healthy muscles in order to improve and facilitate muscle performance
2. Mode 2 (TENS): Relief of pain associated with sore and aching muscles due to strain from exercise or work activities
3. Mode 3 (EMS and TENS): Combination of the benefit from both EMS and TENS of Mode 1 & 2
4. Mode 4 (EMS): Visible muscle contraction saw for improving muscle performance
5. Mode 5 (TENS): Deeper relief of pain associated with sore and aching muscles due to strain from exercise or work activities
6. Mode 6 (EMS and TENS): Combination of the benefit from both EMS and TENS of Mode 4 & 5

MATT MAINTENANCE

- Before applying the unit to the sole feet, make sure to clean the desired area beforehand. To clean skin, simply wet a towel or napkin with rubbing alcohol and gently wipe the area.
- Do not apply lotions or oils before using the device as this will prevent the matt from sticking to the skin.
- The lifespan of the matt will vary depending on the use frequency and care taken. Keeping matt clean and stored properly on the provided plastic sheet will extend their lifespan.
- Keep fluids away from the matt for a longer lifespan. Clean it by swiping the dry tissue or napkin.

TROUBLESHOOTING

- Stimulation is weak or non-existent: Be sure feet are clean and matt are firmly attached.
- Device does not turn on: Recharge for 20 minutes and pin device to the mat properly.

MASSFORCE branded items purchased from massforce.co include a 2-year limited replacement warranty. For Customer Service, please contact us via email hello@massforce.co

The Massforce® name and logo are registered trademarks.
MASSFORCE©. All Rights Reserved.