MASSNECK™ PRO

Modèle LG-A3

MASSFORCE

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MASSFORCE

INSTRUCTION MANUAL

Read the manual before operating this product.
Use only as directed.

WARNINGS

- Do not use this product when there is injury, redness, inflammation, congestion, or other suitable diseases in the neck.
- This product has a heating function. For users who are prone to heat allergies or insensitive to heat sources, please use the heating function with caution.
- If you feel unwell during or after use, please stop using this product.
- This product is a personal care product and cannot replace professional treatment or health care. Please do not use it for diagnosis or other commercial purposes.
- Do not use it with other medical equipment or more than 2 sets of this product at the same time.
- Do not use the same site for more than 30 minutes, otherwise, it may cause discomfort.
- Do not use the heating function for more than 30 minutes, otherwise, it may cause low-temperature skin burns.
- Do not use the product after accidentally falling into water or liquid.

PRECAUTIONS

- This product contains lithium batteries, please do not disassemble, modify, immerse the product in water or throw fire into it by yourself.
- In the event of an abnormal situation that cannot work, please cut off the power immediately and remove it, check the basic troubleshooting or consult customer service.
- When the product is working, do not touch the electrode pads with metal accessories such as belts and necklaces.
- Do not charge, use and place in high-temperature places such as near the stove or under the scorching sun.
- This product is not waterproof, please do not wash and pay attention to waterproof.
- This product can only be decomposed when the product is discarded. If you need to repair it, please consult customer service in time.
- Do not leave the product in front of unattended children.
- Do not allow children or persons with limited ability to use this product is not suitable for any sick or incapacitated person.
- If you have one of the following conditions, please consult a doctor before use:
 - A. Those are pregnant, sick or in poor health,
 - B. Patients who are undergoing treatment and feel abnormalities,
 - C. Patients with high fever and infectious diseases,
 - D. Patients with hypertension, heart disease, cerebral thrombosis, arteriosclerosis or malignant tumors,
 - E. People with skin sensation disorder or abnormal skin sensitivity,
 - F. Those in the postoperative recovery stage,
 - G. Those who are uncertain about their health.

ADVERSE EFFECTS

- Skin irritation and burns beneath the stimulation electrodes applied to your skin can occur.
- On very rare occasions, first-time users of EMS have reported feeling light-headed or faint. We recommend that you use the product while seated until you become accustomed to the sensation.


You should stop using the device and consult with your physician if you experience adverse reactions from the device.

SPECIFICATION

Product model: LG-A3
Product material: ABS, PC, silicone
Lithium battery: 450mAh
Standby current: $100\mu\text{A}$
Input standard: 5V
Maximum rating: 5W
Host size: 105 * 135 * 35mm

***Note:**

This product is not suitable for medical equipment.
Does not have a therapeutic effect,
children can not touch.

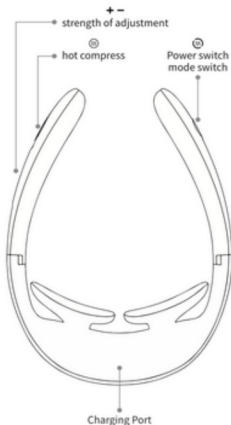INTRODUCTION

MASSNECK™ PRO is the folding portable neck messenger with the material of stainless steel electrode sheets and PI heating sheets, micro-computer chip control technology, and helps relax the neck muscles and relieve the neck fatigue through micro-current stimulation and hot compress. A hot compress can effectively raise the temperature of the environment around the neck, promote blood circulation, and cooperate with the stimulation of microcurrent to drive the neck muscles to exercise together, giving people a full range of comfortable experiences.

MASSNECK™ PRO can be folded to make it easier to carry.

MAINTENANCE

When cleaning and maintaining this product, please turn-off the host and use a wet paper towel to clean the electrode pads. When it is particularly dirty, wipe and disinfect with a little medical alcohol (75% alcohol concentration).

Please use a neutral detergent and soft cloth for daily maintenance. Do not use corrosive detergents. Do not soil or rinse in water. During the cleaning process, water or other liquids should be prevented from entering the body to avoid short circuit and electric shock hazard during use.

USING MASSNECK™

1. Before using this product, please put a hydrogel sheet on it when is turned off, or use a wet towel, wet paper towel or clean water to wipe and moisten your neck massaged area.
2. In order to make the electrode pads fully contact the massage site, please make sure the electrode pads fit the skin tightly when wearing.
3. Do not touch the patch with your hands when it is turned on.

Starting up and shut down.

Press and hold the host button of the product for two seconds, the host will turn on (hear a "di" sound, the blue light is on), enter the default mode, the default gear (0 gear) in the power-on state, press and hold the host button for two seconds, you will hear "di" two beeps, the product shuts down.

Mode selection (currently 5 modes)

After the product host is turned on, short press the host (M) as the mode key, and short press the mode key to switch modes. There are 5 modes in total. Every time you press it, you hear a beep, the mode is incremented once, and the five modes can be switched in cycles (Note: M1 combination mode, M2 massage mode, M3 tapping mode, M4 Gua Sha mode, M5 acupuncture mode). Each time you switch a mode, the velocity returns to the default 0 level.

STRENGTH SELECTION

1. In the power-on state, the "+" on the main unit of the product is the strength strengthening button. Each time you press it, you will hear a beep, the strength will increase by one level, and the strength will be divided into 15 levels (0, 1, 2, ... 15). Press the "+" gear again when the intensity is 15 to make a double beep, the intensity remains unchanged, and the 15 gear is maintained.
2. The "--" on the main unit of the product is the force reduction button. Each time you press it, you hear a beep the intensity decreases by one level, and the force is divided into 15 levels (15, 14, 13, ...0) and then press at 0 in the "--" gear, you will hear two beeps of "Beep", the intensity remains unchanged, and the 0 gear is maintained.

Heating options


1. Short press the (M) button on the main unit of the product to turn on the heating button, press the heating button "di" to turn on the heating function of the main unit, and gradually reach a constant temperature of 42 within 30 seconds.
2. Short press the (M) button on the main unit of the product to turn off the heating button, press the button to turn off the heating "didi" twice to stop the heating function.

Timing function

1. The product host is scheduled for 15 minutes. When the product host is turned on, it will automatically shut down after 15 minutes. It is recommended to use it for no more than 30 minutes a day.
2. When the product host is turned on, the product will automatically shut down within 3 minutes when the product is not worn and used normally.

Low voltage alarm charging function

Low voltage alarm: when the product is out of power, the blue indicator light flashes quickly (frequency 3 times per second), you can plug one end of the charging cable into the product charging port, and the other end into the mobile phone charger with Android head, the red light flashes while charging. The red light is always on when it is full. The product cannot be used while charging.

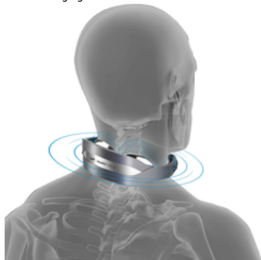TROUBLESHOOTING

If you encounter problems when using the product, please refer to the following tips and solutions.

If the product is still in a faulty state or cannot be used, please turn off the power.

Users are not allowed to disassemble or repair the product by themselves. The product damage caused by disassembly or repair by themselves will invalidate the product warranty, and the company will not bear any consequences caused by this situation.

Unable to start - The internal charging battery is low. Please confirm that the key of the host product is pressed.

Unable to charge or charging a light display - Please check whether the USB cable is connected with the host and the power adapter of the mobile phone.

It shows normal but no sensation - Is the electrode piece in direct contact with the skin? Adjust the position to make it close to the skin. Do you have hair or clothing accessories between the electrodes and the skin? Remove or remove items. Is the force gear 0? Adjust and increase massage strength.

There is a tingling pain - Is the electrode piece close to the neck skin? Adjust the position to make it close to the skin. Is the massage time too long? It is recommended to press 15 minutes once, twice a day. Is the massage too strong, whether the skin is allergic to the electrode piece or not. Patients with a severe allergy should go to the hospital for desensitization treatment first. Cold weather and dry skin? Use wet towel to wet the skin. Is there any wound on the surface of the skin? Do not place the electrode on the wound for massage.

The stimulation is very weak - Is the massage electrode piece close to the neck skin? It is necessary to adjust the position to make it close to the skin and adjust the appropriate force. Is there dirt or oil on the electrode? Clean the electricity Electrode surface. Is the massage strength too big? Can adjust the small strength. Massage gear strength is too small, adjust to increase massage strength.

Feel free to get in touch with us at hello@massforce.co or visit our website massforce.co

Massforce branded items purchased from massforce.co include a 1-year limited replacement warranty. For Customer Service, please contact us via email hello@massforce.co.

DO NOT DISPOSE OF BATTERIES IN FIRE. ALWAYS DISPOSE OF BATTERIES PER LOCAL & FEDERAL GUIDELINES. The Massforce® name and logo are registered trademarks.

©Massforce. All Rights Reserved.

Read Instructions
Before Operating