MASSFORCE

MASSGUN MINI

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Important Safety Advice and Warnings

1. Not for minors.
2. Do not use this product and consult your doctor if you have a muscles injury.
3. Do not use it on a wet surface or unclean skin. It can be used across clothes. Gently press and move on the skin for about 60 seconds.
4. It should not be used above the neck, on bones, or in places with limited muscle coverage. Use this product only on non-painful or unpleasant soft tissues of the body.
5. Do not use if bruising or other signs of pain or discomfort appear.
6. Do not place any object in the product's ventilation.
7. It is forbidden to overflow the product into water, and water must not enter the product's ventilation hole.
8. It is forbidden to put the product into the fire or near a fire source.
9. Do not drop this product intentionally.
10. Do not disassemble the product in any way.
11. Do not put fingers, hair, or other parts of the body close to the axis and back of the massage head for the possibility of getting hurt.

Charging Instructions

1. Connect the USB cord to the product's charging port. Both cables included can be used for charging the device.
2. Connect the USB port to a **5V - 2.4A** min power source (USB charger, computer port...)

Make sure of the power source. If the voltage is different, the device may not be charged and can be damaged.

3. When the charging process begins, a red LED will light up.
- * The product cannot be used during the charging process.
4. The LED will turn green when the battery is fully charged.

Color Indicator for Charging

The color of the light indicates the charging status

Red		No battery left. Need to recharge
Green		Fully charged
Green		Battery capacity ≤10%.
Green		11% - 30% battery capacity
Green		31% - 60% battery capacity
Green		61% - 100% battery capacity

4 speed modes

Higher speed make a stronger massage.

Speed Indicators

-  First speed
-  Second speed
-  Third speed
-  Fourth speed
-  Low battery - device needs to be charged

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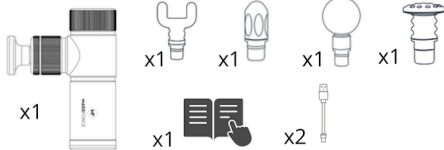
Legs / Shoulder blades

4

1 - "U" shaped massage head: For cervicals and neck.


2 - Cylindrical massage head: For deep tissues and muscles and trigger points (knots).

3 - Round massage head: For large muscle groups such as the back of the thigh and calf, the triceps

4 - Flat massage head: To relax any muscle of the body.

5 - 2 Charging cables USB-C to USB-C / USB to microUSB

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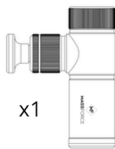
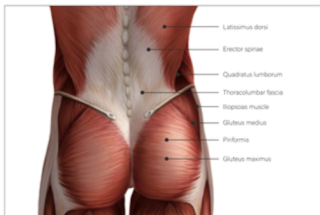
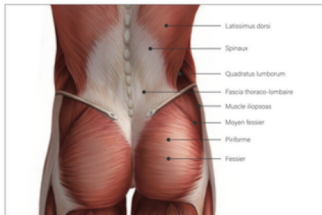
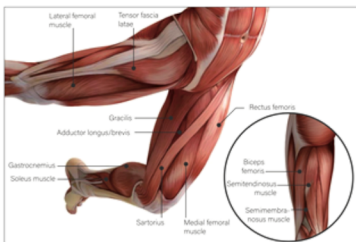
Tips: Massage these body muscles from bottom to top.

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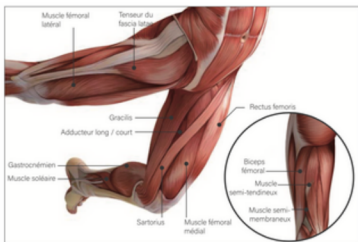Tips: Massage these body muscles from middle to side, from top to bottom.

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

Tips: Massage these body muscles from top to bottom.

ARM

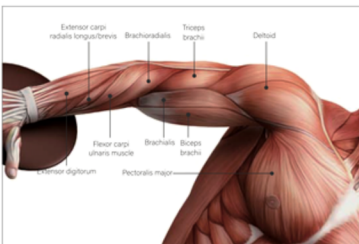Muscle groups	Duration/times of massage
Deltoid	12s-20sx3 times
Biceps brachii	15s-20sx3 times
Triceps brachii	15s-20sx3 times
Brachialis	15s-20sx3 times
Brachioradialis	15s-20sx3 times
Flexor carpi ulnaris muscle	10s-15sx3 times
Extensor carpi radialis longus	5s-10sx2 times
Extensor carpi radialis brevis	5s-10sx2 times
Extensor digitorum	5s-10sx2 times
Pectoralis major	12s-20sx3 times

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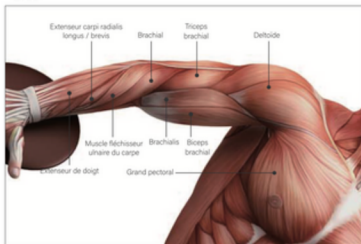USER TIPS

- Move the MassGun Mini for around 1 minute over the desired muscle areas.
- Apply as much pressure as necessary. You should never push too hard since the vibration impulse does all of the work.
- Use the massager for 15-20 minutes each time, up to three times every day.
- Use the massage head in the same area for no more than 3 minutes.
- Excessive usage might cause tissue damage.

CLEANING AND CARE

- Do not immerse the device in water, keep it away from liquids or heat sources, and do not let water to enter the product's vents. Ensure the vents are free of dust and debris.
- Do not drop or improperly use the product, do not remove the screws or attempt to disassemble, do not change or modify the product in any method. Do not operate continuously for more than 1 hour, and let it rest for 30 minutes before reusing it.
- Carefully check the device and battery before use.
- Unplug the cable from the device after charging or before use.
- Do not operate or charge the product unattended.

IMPORTANT: If there is a problem, do not attempt to repair it yourself. There are no warranty claims in such cases.

STORAGE METHOD

Ideally, you should store the Massforce Massgun in its original packaging and keep it in a clean, dry place.

Do not expose to direct sunlight or high temperatures.

If not used for a long period of time, cover it with a dust cloth and remove the power plug from the socket.

Recycling

This symbol indicates that this product may not be disposed of with household waste in any EU country. In order to avoid possible damage to health and the environment, it should be recycled for the reuse of its constituent materials. Recycle the used equipment, follow the recycling process or contact the original dealer.

For the protection of the environment, do not dispose of this product with general household waste at the end of its useful life.

Place the product in the government-designated recycling station, or return it to the manufacturer for service. The service center personnel will handle it in an environmentally friendly manner.

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FABRIQUÉ EN
CHINE

TROUBLESHOOTING

If you encounter problems when using the product, please refer to the following tips and solutions:

1. Please turn off the power if the product is still malfunctioning or cannot be used.
2. Product damage caused by disassembly or repair will invalidate the product warranty, and the company will not be held responsible for any damages resulting from this situation.
3. Unable to start - The battery has low power. Please recharge it.
4. No charging or all of the indication lights are not lighting when charging - the product has a low battery since it has not been used for a long time; fully charge the product.
5. When you notice an insufficient temperature when in warm mode - low battery, fully charge the product.
6. When air pressure mode and warm mode unable to start - Due to a power supply issue, fully charge the product.
7. How to Assess a Charging Situation - When charging, the indicator light remains on. The indicator goes off when the product is completely charged.

SATISFACTION GUARANTEE

If you are dissatisfied with this product for any reason, please CONTACT US at hello@massforce.co and we will make every effort to make it right.

MASSFORCE branded items purchased from massforce.co includes a 1-year limited replacement warranty. For Customer Service, please contact us via email at hello@massforce.co

The MASSFORCE® name and logo are registered trademarks.

©MASSFORCE. All Rights Reserved.

Imported by OUI SMART

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75001 PARIS, FRANCE

www.massforce.co / www.massforce.fr

RoHS

Made in China